

重组蛋白

Catalog Number:ORSMP5406AMu

Recombinant *Macaca fascicularis* Angiotensin-converting
enzyme (ACE2), partial (Active)

Organism Species: *Macaca fascicularis* (Crab-eating
macaque) (*Cynomolgus monkey*)

使用说明手册V1.0

网站: <https://www.aoruisai.com> 手机号: 13061929181(微信)
<https://www.aoruicell.com> 13061929770(微信)

邮箱: sale4@aoruisai.com
sale5@aoruisai.com

电话: 021-57896695
021-57875552

发票/服务/监督投诉电话: 13122716991



企业微信

企业微信

【Product Information】

Gene Names: ACE2

Uniprot:A0A2K5X283

Purity: Greater than 93% as determined by SDS-PAGE.

Endotoxin: Less than 1.0 EU/ug as determined by LAL method.

Source: Mammalian cell

Protein Length: Partial

Expression Region: 18-741aa

Sequence:

QSTIEEQAKTFLDKFNHEAEDLFYQSSLASWNYNTNITEENVQNMNAGEKWSAFLKEQS
TLAQMYPLQEIQNLTVKLQLQALQQNGSSVLSEDKSKRLNTILNTMSTIYSTGKVCNPNN
PQECLLLDPGLNEIMEKSLDYNERLWAWEGWRSEVGKQLRPLYEEYVVLKNEMARANHYK
DYGDYWRGNYEVNGVDGYDYNRDQLIEDVERTFEEIKPLYEHLHAYVRAKLMNAYPSYIS
PTGCLPAHLLGDMWGRFWTNLYSLTVPFQKPNIDVTDAMVNAWNAQRIFKEAEKFFVS
VGLPNMTQGFWENSMLTDPGNVQKVCHPTAWDLGKGDFRIIMCTKVTMDDFLTAHHEMG
HIQYDMAYAAQPFLLRNGANEGFHEAVGEIMSLSAATPKHLKSIGLLSPDFQEDNETEIN
FLLKQALTIVGTLPTFTYMLEKWRWMVFKGEIPKDQWMKKWEMKREIVGVVEPVPHDETY
CDPASLFHVSNDYSFIRYYTRTLYQFQFQEQALCQAAKHEGPLHKCDISNSTEAGQKLLNM
LKLKSEPWTLALENVGAKNMNVRPLLNYFEPLFTWLKDQNKNSFVGWSTDWSPYADQS
IKVRISLKSALGDKAYEWNENMYLFRSSVAYAMRTYFLEIKHQITILFGEEDVRVADLKP
RISFNFYVTAPKNVSDIIPRTEVEEAIRISRSRINDAFRLNDNSLEFLGIQTTLAPPYQS
PVT

Tag Info: C-terminal hFc-tagged

Form:Lyophilized powder

Buffer formulation:Lyophilized from a 0.2 μ m filtered PBS,
6% Trehalose, pH 7.4

Applications:-

Predicted Molecular Mass:-

Accurate Molecular Mass:112.7 kDa

Predicted isoelectric point:-

网站: <https://www.aoruisai.com> 手机号: 13061929181(微信)
<https://www.aoruicell.com> 13061929770(微信)

邮箱: sale4@aoruisai.com
sale5@aoruisai.com

电话: 021-57896695
021-57875552

发票/服务/监督投诉电话: 13122716991




企业微信

企业微信

【Storage and Stability】

Storage: Store at 2-8°C for one month, the shelf life of liquid form is 6 months at -20° C/-80° C, The shelf life of lyophilized form is 12 months at -20° C/-80° C. Avoid repeated freeze/thaw cycles.

【Identification】

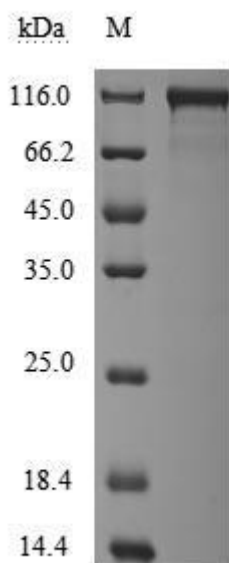

【Important Note】

This product is for scientific research purposes only and should not be used for clinical diagnosis or other procedures, otherwise the consequences will be borne by oneself

网站: <https://www.aoruisai.com> 手机号: 13061929181(微信)
<https://www.aoruicell.com> 13061929770(微信)

邮箱: sale4@aoruisai.com
sale5@aoruisai.com

电话: 021-57896695
021-57875552

发票/服务/监督投诉电话: 13122716991





网站: <https://www.aoruisai.com> 手机号: 13061929181(微信)
<https://www.aoruicell.com> 13061929770(微信)

邮箱: sale4@aoruisai.com
sale5@aoruisai.com

电话: 021-57896695
021-57875552

发票/服务/监督投诉电话: 13122716991



企业微信

企业微信